


Chromatogram : T-01_channel1

System : 450-GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

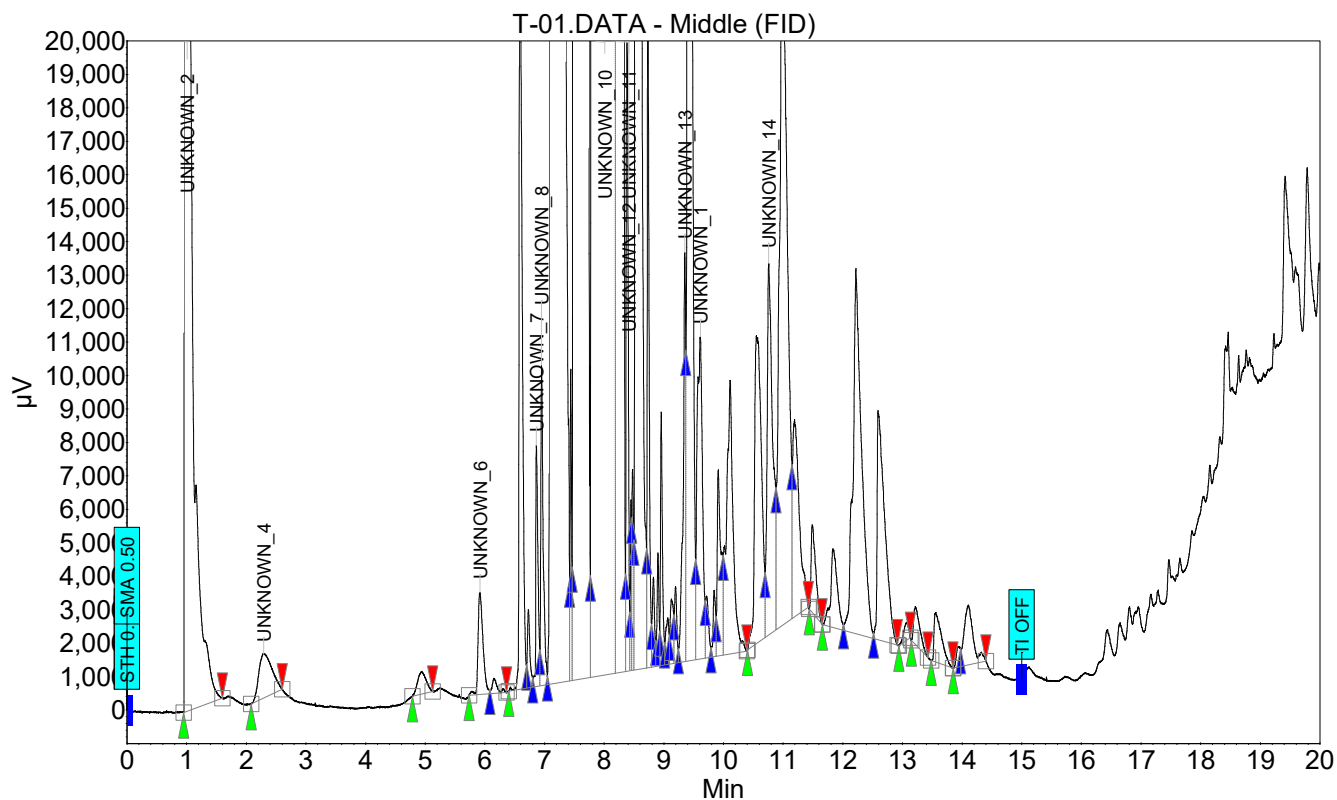

CHROMATOGRAM METHOD REPORT :

Control method :

Scion/Bruker/Varian 400-GC Series (Exclu

Autosampler (8410)

Autosampler	ENABLED
First injector used	Position 2
Syringe volume	10 uL
Advance tray	NO
Clean between injections	NO
Use injection delay	NO
Delay between injections	1.0 min
Injection mode	Std (Split/Splitless)
Sample penetration depth	90 %
Solvent penetration depth	90 %

Abort clean

Chromatogram : T-01_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Vial	I
Volume	5.0 uL
Strokes	1
Drawup speed	5.0 uL/s

Clean mode

Pre-injection solvent flushes	3
Pre-injection sample flushes	0
Post-injection solvent flushes	5
Clean solvent source	I

Injectors

Front (PWOC)

EFC 23 - Flow program	Rate (mL/min/min)	Step (mL/min)	Time (min)
	Initial	10.0	0.00
		Total time	0.00

EFC 23

Enabled OFF

Heat-only zone 1

Heater	OFF
Setpoint	50 °C

Middle (S/SL)

Chromatogram : T-01_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Split event table	Time (min)	Split state	Split ratio
	Initial	ON	25

EFC 25

Enabled ON

Heat-only zone 2

Heater	ON
Setpoint	230 °C

Column oven

Stabilization time 0.50 min

Column oven zone

Heater ON

Temperature program	Rate (°C/min)	Step (°C)	Time (min)
	Initial	100	4.00
	25.0	200	8.00
	20.0	250	5.00
		Total time	23.50

Valve oven

Rear (Small valve oven)

Heat-only zone 3

Heater	OFF
Setpoint	50 °C

Columns

Front

Chromatogram : T-01_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Carrier gas	Helium
Length	15.00 m
Inside diameter	200 um
Type	WCOT

Middle

Carrier gas	Helium
Length	30.00 m
Inside diameter	320 um
Type	PLOT

Pressure mode	Constant flow
Column flow	2.20 mL/min
Pressure pulse	DISABLED

Detectors

Heat-only zone 4

Heater	OFF
Setpoint	50 °C

Front (TCD)

Electronics	OFF
Time constant	Fast
Data rate	10 Hz

Filament temp. limit	390 °C
Filament temperature	50 °C

TCD event table	Time (min)	Range	Autozero	Polarity
	Initial	0.05	YES	POSITIVE

Chromatogram : T-01_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

EFC 14

Enabled	OFF
Make-up (Helium) flow	25 mL/min
Reference (Helium) flow	30.0 mL/min

Heat-only zone 5

Heater	ON
Setpoint	250 °C

Middle (FID)

Electronics	ON
Time constant	Fast
Data rate	10 Hz

FID event table	Time (min)	Range	Autozero
	Initial	12	YES

EFC 11

Enabled	ON
Make-up (Nitrogen) flow	28 mL/min
Combustion (H2) flow	30.0 mL/min
Combustion (Air) flow	300.0 mL/min

Ouput ports

Front

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Front	0

Middle

Chromatogram : T-01_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

ENABLED	NO
---------	----

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Middle	0

Rear

ENABLED	NO
---------	----

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Rear	0

Valves

Valve event table	Time (min)	Gas Sample Valve
	Initial	Fill

Miscellaneous

Start automatically when ready	NO
--------------------------------	----

Acquisition :
Run Name : T-01
Run Id. : -1
Run Time : 24.00

Integration method :
Reduce Noise : Yes
Spike Parameter : 1
Use Relative Threshold ? : Yes
Peak saturation level : min = 0 max = 0

Integration events :

Active	Time	Event	ON	Value
Yes	0.00	Set Peak Width		0.5
Yes	0.00	Set Threshold		0.05
Yes	0.00	Set Minimal Area		0.5
Yes	15.00	Turn Integration	OFF	

Manual actions :

Chromatogram : T-01_channel1

System : 450 GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Manual integration action count : 28

Action #1 : 20/10/2025 03:53:32 p. m. : ADD peak between t=4.787 and t=5.125
Action #2 : 20/10/2025 03:53:50 p. m. : MOVE stop peak # 3 from t=6.090 to t=6.361
Action #3 : 20/10/2025 03:53:54 p. m. : SPLIT peak # 3 at 6.08
Action #4 : 20/10/2025 03:53:56 p. m. : SPLIT peak # 4 at 6.70
Action #5 : 20/10/2025 03:54:32 p. m. : SPLIT peak # 5 at 6.92
Action #6 : 20/10/2025 03:54:51 p. m. : SPLIT peak # 7 at 7.46
Action #7 : 20/10/2025 03:55:11 p. m. : MOVE valley peak # 7 from t=7.457 to t=7.461
Action #8 : 20/10/2025 03:55:27 p. m. : SPLIT peak # 9 at 8.36
Action #9 : 20/10/2025 03:55:36 p. m. : SPLIT peak # 10 at 8.46
Action #10 : 20/10/2025 03:55:40 p. m. : SPLIT peak # 34 at 8.50
Action #11 : 20/10/2025 03:55:42 p. m. : SPLIT peak # 11 at 8.79
Action #12 : 20/10/2025 03:55:45 p. m. : SPLIT peak # 12 at 8.93
Action #13 : 20/10/2025 03:55:47 p. m. : SPLIT peak # 37 at 9.01
Action #14 : 20/10/2025 03:55:49 p. m. : SPLIT peak # 13 at 9.17
Action #15 : 20/10/2025 03:55:51 p. m. : SPLIT peak # 14 at 9.37
Action #16 : 20/10/2025 03:55:57 p. m. : SPLIT peak # 16 at 9.88
Action #17 : 20/10/2025 03:56:09 p. m. : SPLIT peak # 15 at 9.70
Action #18 : 20/10/2025 03:56:59 p. m. : MOVE stop peak # 20 from t=11.233 to t=11.386
Action #19 : 20/10/2025 03:57:19 p. m. : MOVE stop baseline # 5 from t=11.233 y=7607.567 to t=11.389 y=3589.169
Action #20 : 20/10/2025 03:57:31 p. m. : MOVE stop baseline # 5 from t=11.389 y=3589.169 to t=11.425 y=3085.183
Action #21 : 20/10/2025 03:57:35 p. m. : MOVE stop peak # 20 from t=11.386 to t=11.423
Action #22 : 20/10/2025 03:57:42 p. m. : SPLIT peak # 20 at 11.15
Action #23 : 20/10/2025 03:58:43 p. m. : MOVE stop baseline # 7 from t=12.932 y=1926.018 to t=13.134 y=2163.025
Action #24 : 20/10/2025 03:58:47 p. m. : MOVE stop peak # 24 from t=12.932 to t=13.132
Action #25 : 20/10/2025 03:58:57 p. m. : MOVE stop baseline # 7 from t=13.134 y=2163.025 to t=12.918 y=1949.214
Action #26 : 20/10/2025 03:59:02 p. m. : MOVE stop peak # 24 from t=13.132 to t=12.916
Action #27 : 20/10/2025 03:59:06 p. m. : ADD peak between t=12.942 and t=13.134
Action #28 : 20/10/2025 03:59:11 p. m. : SPLIT peak # 27 at 13.98

Peak Identification table :

Peak Name	RT [min]	Abs.Window[min]	Window %	Ref?	Mode
UNKNOWN_2	1.14	0.20	0.00		Nearest
heptano	1.48	0.20	0.00		Nearest
UNKNOWN_4	2.11	0.20	0.00		Nearest
UNKNOWN_5	4.10	0.20	0.00		Nearest
UNKNOWN_6	5.86	0.20	0.00		Nearest
UNKNOWN_7	6.93	0.20	0.00		Nearest
UNKNOWN_8	6.96	0.20	0.00		Nearest
UNKNOWN_9	7.71	0.20	0.00		Nearest
UNKNOWN_10	7.84	0.20	0.00		Nearest
UNKNOWN_11	8.04	0.20	0.00		Nearest
UNKNOWN_12	8.41	0.20	0.00		Nearest
UNKNOWN_13	9.32	0.20	0.00		Nearest
UNKNOWN_1	9.54	0.20	0.00		Nearest
UNKNOWN_14	10.71	0.20	0.00		Nearest

Resolve with references : No

Group Identification table :

Group Name	Group type	Parameters
------------	------------	------------

Chromatogram : T-01_channel1

System : 450-GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Peak results :

Index	Name	Time [Min]	Quantity [% Area]	Height [µV]	Area [µV.Min]	Area % [%]
1	UNKNOWN 2	0.97	79.76	571172335.4	4196392.2	79.759
2	UNKNOWN 4	2.29	0.01	1304.0	336.3	0.006
28	UNKNOWN	4.94	0.00	660.5	111.6	0.002
3	UNKNOWN 6	5.92	0.01	3040.4	292.9	0.006
29	UNKNOWN	6.15	0.00	439.9	43.9	0.001
4	UNKNOWN	6.59	0.03	24815.0	1426.9	0.027
30	UNKNOWN	6.73	0.00	2339.7	116.8	0.002
5	UNKNOWN 7	6.86	0.01	7199.4	307.4	0.006
31	UNKNOWN 8	6.95	0.01	9582.5	450.1	0.009
6	UNKNOWN	7.09	1.61	742009.7	84661.9	1.609
7	UNKNOWN	7.44	0.00	9309.6	246.7	0.005
32	UNKNOWN	7.48	1.49	598967.7	78506.2	1.492
8	UNKNOWN 10	7.87	15.56	3666476.7	818799.0	15.563
9	UNKNOWN 11	8.20	0.85	666529.7	44624.7	0.848
33	UNKNOWN 12	8.38	0.02	22498.9	846.6	0.016
10	UNKNOWN	8.45	0.00	5106.5	111.7	0.002
34	UNKNOWN	8.48	0.00	5983.1	195.5	0.004
35	UNKNOWN	8.53	0.24	171940.1	12723.7	0.242
11	UNKNOWN	8.73	0.01	21576.4	749.0	0.014
36	UNKNOWN	8.83	0.00	2677.5	118.9	0.002
12	UNKNOWN	8.90	0.00	3343.6	126.5	0.002
37	UNKNOWN	8.96	0.01	7574.7	282.7	0.005
38	UNKNOWN	9.07	0.00	1378.4	86.3	0.002
13	UNKNOWN	9.13	0.00	1938.4	117.8	0.002
39	UNKNOWN	9.20	0.00	2286.0	111.1	0.002
14	UNKNOWN 13	9.35	0.01	12226.2	589.9	0.011
40	UNKNOWN	9.39	0.11	80098.3	5714.8	0.109
15	UNKNOWN 1	9.61	0.02	9636.6	997.9	0.019
42	UNKNOWN	9.72	0.00	1846.2	108.3	0.002
16	UNKNOWN	9.84	0.00	1976.4	96.1	0.002
41	UNKNOWN	9.91	0.01	5567.5	415.8	0.008
17	UNKNOWN	10.11	0.02	8180.3	1054.3	0.020
18	UNKNOWN	10.55	0.02	9204.7	1230.9	0.023
19	UNKNOWN 14	10.76	0.02	11107.9	1197.9	0.023
20	UNKNOWN	10.98	0.06	21933.3	3247.2	0.062
43	UNKNOWN	11.19	0.01	5891.3	768.0	0.015
21	UNKNOWN	11.49	0.00	2627.3	231.3	0.004
22	UNKNOWN	11.83	0.01	2340.9	311.6	0.006
23	UNKNOWN	12.22	0.03	10918.6	1829.1	0.035
24	UNKNOWN	12.59	0.02	6844.2	1010.6	0.019
44	UNKNOWN	13.07	0.00	530.0	50.2	0.001
25	UNKNOWN	13.22	0.00	1142.2	137.2	0.003
26	UNKNOWN	13.56	0.00	1484.5	214.9	0.004
27	UNKNOWN	13.94	0.00	599.2	47.0	0.001
45	UNKNOWN	14.11	0.01	1777.1	311.9	0.006
Total			100.00	577347246.4	5261351.3	100.000

CHROMATOGRAM METHOD REPORT :

Control method :

Scion/Bruker/Varian 400-GC Series (Exclu

Chromatogram : T-01_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Autosampler (8410)

Autosampler	ENABLED
First injector used	Position 2
Syringe volume	10 uL
Advance tray	NO
Clean between injections	NO
Use injection delay	NO
Delay between injections	1.0 min
Injection mode	Std (Split/Splitless)
Sample penetration depth	90 %
Solvent penetration depth	90 %

Abort clean

Vial	I
Volume	5.0 uL
Strokes	1
Drawup speed	5.0 uL/s

Clean mode

Pre-injection solvent flushes	3
Pre-injection sample flushes	0
Post-injection solvent flushes	5
Clean solvent source	I

Injectors

Front (PWOC)

EFC 23 - Flow program	Rate (mL/min/min)	Step (mL/min)	Time (min)
	Initial	10.0	0.00
		Total time	0.00

EFC 23

Chromatogram : T-01_channel1

System : 450 GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Enabled	OFF
---------	-----

Heat-only zone 1

Heater	OFF
Setpoint	50 °C

Middle (S/SL)

Split event table	Time (min)	Split state	Split ratio
	Initial	ON	25

EFC 25

Enabled	ON
---------	----

Heat-only zone 2

Heater	ON
Setpoint	230 °C

Column oven

Stabilization time	0.50 min
--------------------	----------

Column oven zone

Heater	ON
--------	----

Temperature program	Rate (°C/min)	Step (°C)	Time (min)
	Initial	100	4.00
	25.0	200	8.00
	20.0	250	5.00
		Total time	23.50

Chromatogram : T-01_channel1

System : 450 GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Valve oven

Rear (Small valve oven)

Heat-only zone 3

Heater	OFF
Setpoint	50 °C

Columns

Front

Carrier gas	Helium
Length	15.00 m
Inside diameter	200 um
Type	WCOT

Middle

Carrier gas	Helium
Length	30.00 m
Inside diameter	320 um
Type	PLOT

Pressure mode	Constant flow
Column flow	2.20 mL/min
Pressure pulse	DISABLED

Detectors

Heat-only zone 4

Heater	OFF
Setpoint	50 °C

Front (TCD)

Electronics	OFF
Time constant	Fast

Chromatogram : T-01_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Electronics	OFF
Data rate	10 Hz

Filament temp. limit	390 °C
Filament temperature	50 °C

TCD event table	Time (min)	Range	Autozero	Polarity
	Initial	0.05	YES	POSITIVE

EFC 14

Enabled	OFF
Make-up (Helium) flow	25 mL/min
Reference (Helium) flow	30.0 mL/min

Heat-only zone 5

Heater	ON
Setpoint	250 °C

Middle (FID)

Electronics	ON
Time constant	Fast
Data rate	10 Hz

FID event table	Time (min)	Range	Autozero
	Initial	12	YES

EFC 11

Chromatogram : T-01_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Enabled	ON
Make-up (Nitrogen) flow	28 mL/min
Combustion (H2) flow	30.0 mL/min
Combustion (Air) flow	300.0 mL/min

Output ports

Front

ENABLED	NO
---------	----

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Front	0

Middle

ENABLED	NO
---------	----

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Middle	0

Rear

ENABLED	NO
---------	----

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Rear	0

Valves

Valve event table	Time (min)	Gas Sample Valve
	Initial	Fill

Miscellaneous

Start automatically when ready	NO
--------------------------------	----

Chromatogram : T-01_channel1

System : 456-GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Acquisition :
System : 456-GC
Project : Aceites
Run Name : T-01
Description
Run Id. : -1
Run Time : 24.00
Vial : 4
Rack : 0
Analysis : Unknown
Injection volume : 1.00

Scale :
RT Min : 0.00
RT Max : 20.00
Y Min : -1,000.00
Y Max : 20,000.00
Divisor factor : 1.00
Multiplier factor : 1.00
Sample mass : 0.00

Run log :

Injection report for Position 2

Injection occurred at 20/10/2025 11:41:23 a. m.

Scion/Bruker/Varian 400-GC Series (Exclu [Scion/Bruker/Varian 400-GC Series (Excluding 430)]
Wrapper 5.0.8.31952
Firmware 32834
Driver 4.1.0.394

Preprocessing :
Blank Subtract :
File:
N.A.

Integration method :
Reduce Noise : Yes
Spike Parameter : 1
Use Relative Threshold ? : Yes
Peak saturation level : min = 0 max = 0

Integration events :

Active	Time	Event	ON	Value
Yes	0.00	Set Peak Width		0.5
Yes	0.00	Set Threshold		0.05
Yes	0.00	Set Minimal Area		0.5
Yes	15.00	Turn Integration	OFF	

Manual actions :
Manual integration action count : 28
Action #1 : 20/10/2025 03:53:32 p. m. : ADD peak between t=4.787 and t=5.125
Action #2 : 20/10/2025 03:53:50 p. m. : MOVE stop peak # 3 from t=6.090 to t=6.361
Action #3 : 20/10/2025 03:53:54 p. m. : SPLIT peak # 3 at 6.08

Chromatogram : T-01_channel1

System : 450 GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Action #4 : 20/10/2025 03:53:56 p. m. : SPLIT peak # 4 at 6.70
Action #5 : 20/10/2025 03:54:32 p. m. : SPLIT peak # 5 at 6.92
Action #6 : 20/10/2025 03:54:51 p. m. : SPLIT peak # 7 at 7.46
Action #7 : 20/10/2025 03:55:11 p. m. : MOVE valley peak # 7 from t=7.457 to t=7.461
Action #8 : 20/10/2025 03:55:27 p. m. : SPLIT peak # 9 at 8.36
Action #9 : 20/10/2025 03:55:36 p. m. : SPLIT peak # 10 at 8.46
Action #10 : 20/10/2025 03:55:40 p. m. : SPLIT peak # 34 at 8.50
Action #11 : 20/10/2025 03:55:42 p. m. : SPLIT peak # 11 at 8.79
Action #12 : 20/10/2025 03:55:45 p. m. : SPLIT peak # 12 at 8.93
Action #13 : 20/10/2025 03:55:47 p. m. : SPLIT peak # 37 at 9.01
Action #14 : 20/10/2025 03:55:49 p. m. : SPLIT peak # 13 at 9.17
Action #15 : 20/10/2025 03:55:51 p. m. : SPLIT peak # 14 at 9.37
Action #16 : 20/10/2025 03:55:57 p. m. : SPLIT peak # 16 at 9.88
Action #17 : 20/10/2025 03:56:09 p. m. : SPLIT peak # 15 at 9.70
Action #18 : 20/10/2025 03:56:59 p. m. : MOVE stop peak # 20 from t=11.233 to t=11.386
Action #19 : 20/10/2025 03:57:19 p. m. : MOVE stop baseline # 5 from t=11.233 y=7607.567 to t=11.389 y=3589.169
Action #20 : 20/10/2025 03:57:31 p. m. : MOVE stop baseline # 5 from t=11.389 y=3589.169 to t=11.425 y=3085.183
Action #21 : 20/10/2025 03:57:35 p. m. : MOVE stop peak # 20 from t=11.386 to t=11.423
Action #22 : 20/10/2025 03:57:42 p. m. : SPLIT peak # 20 at 11.15
Action #23 : 20/10/2025 03:58:43 p. m. : MOVE stop baseline # 7 from t=12.932 y=1926.018 to t=13.134 y=2163.025
Action #24 : 20/10/2025 03:58:47 p. m. : MOVE stop peak # 24 from t=12.932 to t=13.132
Action #25 : 20/10/2025 03:58:57 p. m. : MOVE stop baseline # 7 from t=13.134 y=2163.025 to t=12.918 y=1949.214
Action #26 : 20/10/2025 03:59:02 p. m. : MOVE stop peak # 24 from t=13.132 to t=12.916
Action #27 : 20/10/2025 03:59:06 p. m. : ADD peak between t=12.942 and t=13.134
Action #28 : 20/10/2025 03:59:11 p. m. : SPLIT peak # 27 at 13.98

Peak Identification table :

Peak Name	RT [min]	Abs.Window[min]	Window %	Ref?	Mode
UNKNOWN_2	1.14	0.20	0.00		Nearest
heptano	1.48	0.20	0.00		Nearest
UNKNOWN_4	2.11	0.20	0.00		Nearest
UNKNOWN_5	4.10	0.20	0.00		Nearest
UNKNOWN_6	5.86	0.20	0.00		Nearest
UNKNOWN_7	6.93	0.20	0.00		Nearest
UNKNOWN_8	6.96	0.20	0.00		Nearest
UNKNOWN_9	7.71	0.20	0.00		Nearest
UNKNOWN_10	7.84	0.20	0.00		Nearest
UNKNOWN_11	8.04	0.20	0.00		Nearest
UNKNOWN_12	8.41	0.20	0.00		Nearest
UNKNOWN_13	9.32	0.20	0.00		Nearest
UNKNOWN_1	9.54	0.20	0.00		Nearest
UNKNOWN_14	10.71	0.20	0.00		Nearest

Resolve with references : No

Group Identification table :

Group Name	Group type	Parameters
------------	------------	------------

Calibration Method :

Method type : Response %

Response : Area

Standard Unit : g/l

Calibration curve file name :

Factors : Curve

Subtract internal standard mass

Response unit : Curve

Chromatogram : T-01_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 20/10/2025 11:41:23 a. m.
Processed : 20/10/2025 03:59:14 p. m.
Printed : 20/10/2025 04:01:24 p. m.

Unknown mode : None

Calibration components :

Component	Level 1	Level 2	Level 3	Level 4	Control Sample
-----------	---------	---------	---------	---------	----------------

Suitability tests :

Export :

Post Run :

Print Method :
File name : : default_standard
Copies : 0

Summary :