


Chromatogram : 6.2_channel1

System : 450-GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

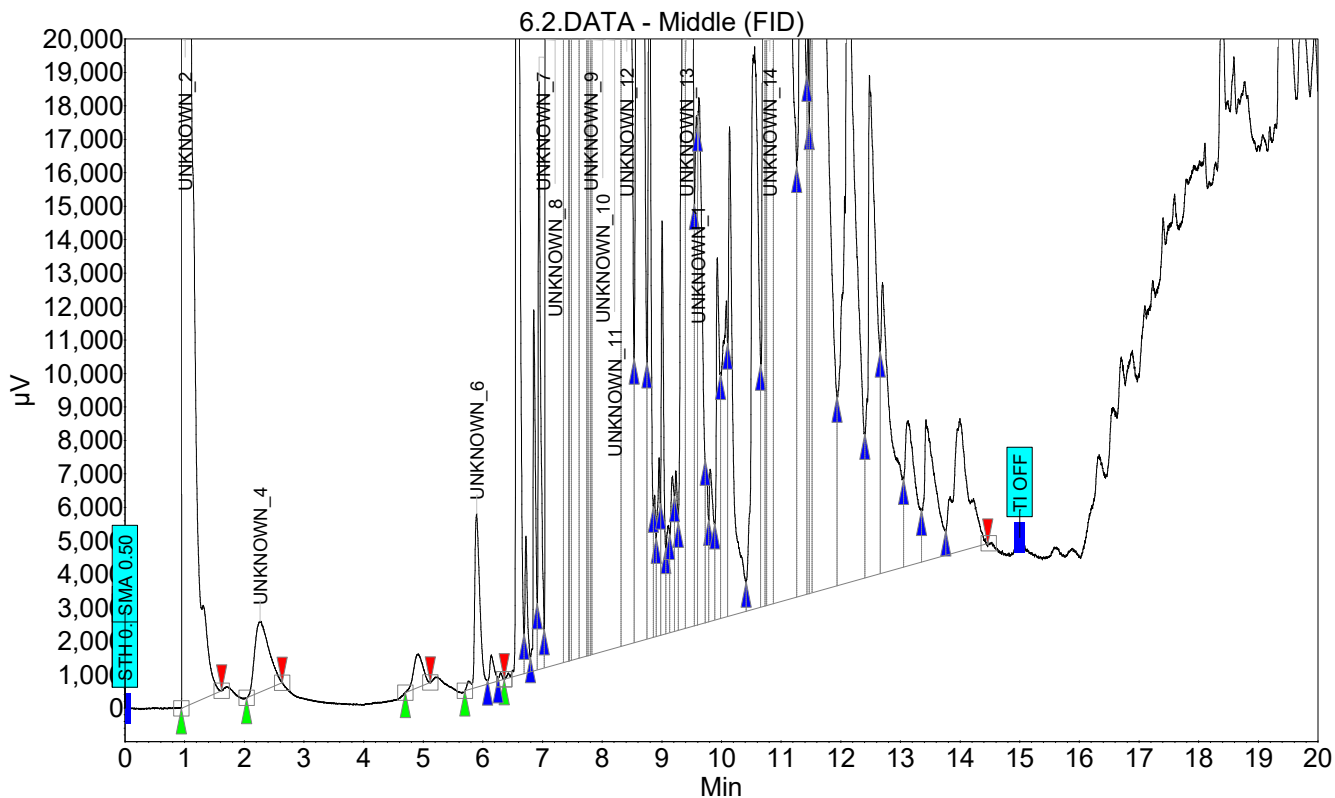

CHROMATOGRAM METHOD REPORT :

Control method :

Scion/Bruker/Varian 400-GC Series (Exclu

Autosampler (8410)

Autosampler	ENABLED
First injector used	Position 2
Syringe volume	10 uL
Advance tray	NO
Clean between injections	NO
Use injection delay	NO
Delay between injections	1.0 min
Injection mode	Std (Split/Splitless)
Sample penetration depth	90 %
Solvent penetration depth	90 %

Abort clean

Chromatogram : 6.2_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

Vial	I
Volume	5.0 uL
Strokes	1
Drawup speed	5.0 uL/s

Clean mode

Pre-injection solvent flushes	3
Pre-injection sample flushes	0
Post-injection solvent flushes	5
Clean solvent source	I

Injectors

Front (PWOC)

EFC 23 - Flow program	Rate (mL/min/min)	Step (mL/min)	Time (min)
	Initial	10.0	0.00
		Total time	0.00

EFC 23

Enabled OFF

Heat-only zone 1

Heater	OFF
Setpoint	50 °C

Middle (S/SL)

Chromatogram : 6.2_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

Split event table	Time (min)	Split state	Split ratio
	Initial	ON	25

EFC 25

Enabled ON

Heat-only zone 2

Heater	ON
Setpoint	230 °C

Column oven

Stabilization time 0.50 min

Column oven zone

Heater ON

Temperature program	Rate (°C/min)	Step (°C)	Time (min)
	Initial	100	4.00
	25.0	200	8.00
	20.0	250	5.00
		Total time	23.50

Valve oven

Rear (Small valve oven)

Heat-only zone 3

Heater	OFF
Setpoint	50 °C

Columns

Front

Chromatogram : 6.2_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

Carrier gas	Helium
Length	15.00 m
Inside diameter	200 um
Type	WCOT

Middle

Carrier gas	Helium
Length	30.00 m
Inside diameter	320 um
Type	PLOT

Pressure mode	Constant flow
Column flow	2.20 mL/min
Pressure pulse	DISABLED

Detectors

Heat-only zone 4

Heater	OFF
Setpoint	50 °C

Front (TCD)

Electronics	OFF
Time constant	Fast
Data rate	10 Hz

Filament temp. limit	390 °C
Filament temperature	50 °C

TCD event table	Time (min)	Range	Autozero	Polarity
	Initial	0.05	YES	POSITIVE

Chromatogram : 6.2_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

EFC 14

Enabled	OFF
Make-up (Helium) flow	25 mL/min
Reference (Helium) flow	30.0 mL/min

Heat-only zone 5

Heater	ON
Setpoint	250 °C

Middle (FID)

Electronics	ON
Time constant	Fast
Data rate	10 Hz

FID event table	Time (min)	Range	Autozero
	Initial	12	YES

EFC 11

Enabled	ON
Make-up (Nitrogen) flow	28 mL/min
Combustion (H2) flow	30.0 mL/min
Combustion (Air) flow	300.0 mL/min

Ouput ports

Front

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Front	0

Middle

Chromatogram : 6.2_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

ENABLED	NO
---------	----

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Middle	0

Rear

ENABLED	NO
---------	----

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Rear	0

Valves

Valve event table	Time (min)	Gas Sample Valve
	Initial	Fill

Miscellaneous

Start automatically when ready	NO
--------------------------------	----

Acquisition :
Run Name : 6.2
Run Id. : -1
Run Time : 24.00

Integration method :
Reduce Noise : Yes
Spike Parameter : 1
Use Relative Threshold ? : Yes
Peak saturation level : min = 0 max = 0

Integration events :

Active	Time	Event	ON	Value
Yes	0.00	Set Peak Width		0.5
Yes	0.00	Set Threshold		0.05
Yes	0.00	Set Minimal Area		0.5
Yes	15.00	Turn Integration	OFF	

Manual actions :

Chromatogram : 6.2_channel1

System : 456 GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

Manual integration action count : 29

Action #1 : 21/11/2025 12:34:24 p. m. : MOVE stop baseline # 5 from t=14.333 y=5443.809 to t=14.478 y=4924.442
Action #2 : 21/11/2025 12:34:24 p. m. : MOVE stop peak # 31 from t=14.333 to t=14.340
Action #3 : 21/11/2025 12:34:24 p. m. : MOVE stop baseline # 5 from t=14.478 y=4924.442 to t=14.478 y=4924.442
Action #4 : 21/11/2025 12:34:24 p. m. : MOVE stop peak # 31 from t=14.340 to t=14.465
Action #5 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 5 at 6.26
Action #6 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 6 at 6.69
Action #7 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 7 at 6.91
Action #8 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 10 at 7.48
Action #9 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 14 at 8.86
Action #10 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 15 at 8.98
Action #11 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 16 at 9.17
Action #12 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 18 at 9.73
Action #13 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 21 at 10.10
Action #14 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 23 at 10.73
Action #15 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 41 at 10.76
Action #16 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 24 at 11.43
Action #17 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 25 at 11.52
Action #18 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 12 at 7.80
Action #19 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 45 at 7.83
Action #20 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 11 at 7.74
Action #21 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 46 at 8.32
Action #22 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 17 at 9.39
Action #23 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 42 at 10.87
Action #24 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 8 at 7.35
Action #25 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 16 at 9.12
Action #26 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 38 at 9.21
Action #27 : 21/11/2025 12:34:24 p. m. : MOVE valley peak # 16 from t=9.123 to t=9.132
Action #28 : 21/11/2025 12:34:24 p. m. : MERGE peak # 52
Action #29 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 18 at 9.60

Peak Identification table :

Peak Name	RT [min]	Abs.Window[min]	Window %	Ref?	Mode
UNKNOWN_2	1.14	0.20	0.00		Nearest
heptano	1.48	0.20	0.00		Nearest
UNKNOWN_4	2.11	0.20	0.00		Nearest
UNKNOWN_5	4.10	0.20	0.00		Nearest
UNKNOWN_6	5.86	0.20	0.00		Nearest
UNKNOWN_7	6.93	0.20	0.00		Nearest
UNKNOWN_8	6.96	0.20	0.00		Nearest
UNKNOWN_9	7.71	0.20	0.00		Nearest
UNKNOWN_10	7.84	0.20	0.00		Nearest
UNKNOWN_11	8.04	0.20	0.00		Nearest
UNKNOWN_12	8.41	0.20	0.00		Nearest
UNKNOWN_13	9.32	0.20	0.00		Nearest
UNKNOWN_1	9.54	0.20	0.00		Nearest
UNKNOWN_14	10.71	0.20	0.00		Nearest

Resolve with references : No

Group Identification table :

Group Name	Group type	Parameters
------------	------------	------------

Chromatogram : 6.2_channel1

System : 450 GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

Peak results :

Index	Name	Time [Min]	Quantity [% Area]	Height [µV]	Area [µV.Min]	Area % [%]
1	UNKNOWN 2	0.96	55.80	394168372.4	2617156.8	55.804
2	UNKNOWN 4	2.26	0.01	2106.3	584.7	0.012
3	UNKNOWN	4.91	0.00	990.6	179.4	0.004
4	UNKNOWN 6	5.89	0.01	5183.9	519.3	0.011
5	UNKNOWN	6.14	0.00	835.7	71.0	0.002
32	UNKNOWN	6.30	0.00	208.6	11.6	0.000
6	UNKNOWN	6.56	0.06	44824.5	2935.1	0.063
33	UNKNOWN	6.72	0.00	4094.7	207.0	0.004
7	UNKNOWN	6.85	0.01	10779.6	492.4	0.010
34	UNKNOWN 7	6.94	0.02	17575.5	880.7	0.019
8	UNKNOWN 8	7.05	3.36	1184025.8	157356.2	3.355
51	UNKNOWN	7.39	0.22	145176.9	10486.1	0.224
9	UNKNOWN	7.43	0.02	75412.7	991.1	0.021
10	UNKNOWN	7.47	0.05	91391.5	2162.7	0.046
35	UNKNOWN	7.50	0.89	511590.4	41814.5	0.892
11	UNKNOWN 9	7.71	0.83	521939.2	38719.5	0.826
47	UNKNOWN	7.75	0.12	568033.7	5748.7	0.123
12	UNKNOWN	7.78	1.67	2701187.5	78350.8	1.671
45	UNKNOWN 10	7.82	1.57	2833290.3	73607.2	1.569
46	UNKNOWN 11	8.07	32.11	4641008.8	1506058.9	32.112
48	UNKNOWN 12	8.33	1.01	574354.4	47470.6	1.012
13	UNKNOWN	8.60	0.51	315964.2	24140.3	0.515
14	UNKNOWN	8.79	0.03	32266.7	1591.9	0.034
36	UNKNOWN	8.88	0.00	4235.0	169.7	0.004
15	UNKNOWN	8.95	0.01	5319.5	315.1	0.007
37	UNKNOWN	9.00	0.01	12367.5	606.3	0.013
16	UNKNOWN	9.11	0.00	3224.0	194.6	0.004
38	UNKNOWN	9.18	0.01	4660.0	324.5	0.007
52	UNKNOWN	9.24	0.01	4768.6	258.7	0.006
17	UNKNOWN 13	9.37	0.09	90657.6	4374.4	0.093
49	UNKNOWN	9.40	0.23	139621.4	10968.0	0.234
18	UNKNOWN 1	9.58	0.02	15240.2	767.1	0.016
53	UNKNOWN	9.62	0.03	15749.9	1413.4	0.030
39	UNKNOWN	9.73	0.01	4851.1	247.6	0.005
19	UNKNOWN	9.81	0.01	4547.9	373.6	0.008
20	UNKNOWN	9.93	0.02	10812.0	711.7	0.015
21	UNKNOWN	10.08	0.02	9443.2	983.6	0.021
40	UNKNOWN	10.13	0.03	14615.5	1418.6	0.030
22	UNKNOWN	10.56	0.05	16794.6	2362.8	0.050
23	UNKNOWN 14	10.71	0.03	31072.6	1592.8	0.034
41	UNKNOWN	10.75	0.02	35247.6	849.6	0.018
42	UNKNOWN	10.85	0.13	66074.5	5941.2	0.127
50	UNKNOWN	10.89	0.39	119601.2	18164.5	0.387
24	UNKNOWN	11.31	0.07	25159.1	3272.6	0.070
43	UNKNOWN	11.45	0.01	18081.6	625.4	0.013
25	UNKNOWN	11.50	0.02	24584.1	888.8	0.019
44	UNKNOWN	11.54	0.20	44099.9	9331.2	0.199
26	UNKNOWN	12.12	0.11	20092.2	4964.6	0.106
27	UNKNOWN	12.48	0.05	14980.6	2427.8	0.052
28	UNKNOWN	12.70	0.04	8684.4	1885.5	0.040
29	UNKNOWN	13.13	0.02	4342.4	889.5	0.019
30	UNKNOWN	13.43	0.02	4221.3	914.3	0.019
31	UNKNOWN	14.00	0.03	3979.1	1177.0	0.025

Chromatogram : 6.2_channel1

System : 450-GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

Index	Name	Time [Min]	Quantity [% Area]	Height [μV]	Area [μV.Min]	Area % [%]
Total			100.00	409227742.6	4689950.7	100.000

CHROMATOGRAM METHOD REPORT :

Control method :

Scion/Bruker/Varian 400-GC Series (Exclu

Autosampler (8410)

Autosampler	ENABLED
First injector used	Position 2
Syringe volume	10 uL
Advance tray	NO
Clean between injections	NO
Use injection delay	NO
Delay between injections	1.0 min
Injection mode	Std (Split/Splitless)
Sample penetration depth	90 %
Solvent penetration depth	90 %

Abort clean

Vial	I
Volume	5.0 uL
Strokes	1
Drawup speed	5.0 uL/s

Clean mode

Pre-injection solvent flushes	3
Pre-injection sample flushes	0
Post-injection solvent flushes	5
Clean solvent source	I

Injectors

Front (PWOC)

Chromatogram : 6.2_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

EFC 23 - Flow program	Rate (mL/min/min)	Step (mL/min)	Time (min)
	Initial	10.0	0.00
		Total time	0.00

EFC 23

Enabled OFF

Heat-only zone 1

Heater	OFF
Setpoint	50 °C

Middle (S/SL)

Split event table	Time (min)	Split state	Split ratio
	Initial	ON	25

EFC 25

Enabled ON

Heat-only zone 2

Heater	ON
Setpoint	230 °C

Column oven

Stabilization time 0.50 min

Column oven zone

Chromatogram : 6.2_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

Heater	ON
--------	----

Temperature program	Rate (°C/min)	Step (°C)	Time (min)
	Initial	100	4.00
	25.0	200	8.00
	20.0	250	5.00
		Total time	23.50

Valve oven

Rear (Small valve oven)

Heat-only zone 3

Heater	OFF
Setpoint	50 °C

Columns

Front

Carrier gas	Helium
Length	15.00 m
Inside diameter	200 um
Type	WCOT

Middle

Carrier gas	Helium
Length	30.00 m
Inside diameter	320 um
Type	PLOT

Pressure mode	Constant flow
Column flow	2.20 mL/min
Pressure pulse	DISABLED

Detectors

Chromatogram : 6.2_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

Heat-only zone 4

Heater	OFF
Setpoint	50 °C

Front (TCD)

Electronics	OFF
Time constant	Fast
Data rate	10 Hz

Filament temp. limit	390 °C
Filament temperature	50 °C

TCD event table	Time (min)	Range	Autozero	Polarity
	Initial	0.05	YES	POSITIVE

EFC 14

Enabled	OFF
Make-up (Helium) flow	25 mL/min
Reference (Helium) flow	30.0 mL/min

Heat-only zone 5

Heater	ON
Setpoint	250 °C

Middle (FID)

Electronics	ON
Time constant	Fast
Data rate	10 Hz

FID event table	Time (min)	Range	Autozero
-----------------	------------	-------	----------

Chromatogram : 6.2_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

FID event table	Time (min)	Range	Autozero
	Initial	12	YES

EFC 11

Enabled	ON
Make-up (Nitrogen) flow	28 mL/min
Combustion (H2) flow	30.0 mL/min
Combustion (Air) flow	300.0 mL/min

Ouput ports

Front

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Front	0

Middle

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Middle	0

Rear

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Rear	0

Valves

Valve event table	Time (min)	Gas Sample Valve
	Initial	Fill

Chromatogram : 6.2_channel1

System : 456-GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

Miscellaneous

Start automatically when ready	NO
--------------------------------	----

Acquisition :
System : 456-GC
Project : Aceites
Run Name : 6.2
Description
Run Id. : -1
Run Time : 24.00
Vial : 6
Rack : 0
Analysis : Unknown
Injection volume : 1.00

Scale :
RT Min : 0.00
RT Max : 20.00
Y Min : -1,000.00
Y Max : 20,000.00
Divisor factor : 1.00
Multiplier factor : 1.00
Sample mass : 0.00

Run log :

Injection report for Position 2

Injection occurred at 07/11/2025 11:34:48 a. m.

Scion/Bruker/Varian 400-GC Series (Exclu [Scion/Bruker/Varian 400-GC Series (Excluding 430)]
Wrapper 5.0.8.31952
Firmware 32834
Driver 4.1.0.394

Preprocessing :
Blank Subtract :
File:
N.A.

Integration method :
Reduce Noise : Yes
Spike Parameter : 1
Use Relative Threshold ? : Yes
Peak saturation level : min = 0 max = 0

Integration events :

Active	Time	Event	ON	Value
--------	------	-------	----	-------

Chromatogram : 6.2_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

Active	Time	Event	ON	Value
Yes	0.00	Set Peak Width		0.5
Yes	0.00	Set Threshold		0.05
Yes	0.00	Set Minimal Area		0.5
Yes	15.00	Turn Integration	OFF	

Manual actions :

Manual integration action count : 29

Action #1 : 21/11/2025 12:34:24 p. m. : MOVE stop baseline # 5 from t=14.333 y=5443.809 to t=14.478 y=4924.442
Action #2 : 21/11/2025 12:34:24 p. m. : MOVE stop peak # 31 from t=14.333 to t=14.340
Action #3 : 21/11/2025 12:34:24 p. m. : MOVE stop baseline # 5 from t=14.478 y=4924.442 to t=14.478 y=4924.442
Action #4 : 21/11/2025 12:34:24 p. m. : MOVE stop peak # 31 from t=14.340 to t=14.465
Action #5 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 5 at 6.26
Action #6 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 6 at 6.69
Action #7 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 7 at 6.91
Action #8 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 10 at 7.48
Action #9 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 14 at 8.86
Action #10 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 15 at 8.98
Action #11 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 16 at 9.17
Action #12 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 18 at 9.73
Action #13 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 21 at 10.10
Action #14 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 23 at 10.73
Action #15 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 41 at 10.76
Action #16 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 24 at 11.43
Action #17 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 25 at 11.52
Action #18 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 12 at 7.80
Action #19 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 45 at 7.83
Action #20 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 11 at 7.74
Action #21 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 46 at 8.32
Action #22 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 17 at 9.39
Action #23 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 42 at 10.87
Action #24 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 8 at 7.35
Action #25 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 16 at 9.12
Action #26 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 38 at 9.21
Action #27 : 21/11/2025 12:34:24 p. m. : MOVE valley peak # 16 from t=9.123 to t=9.132
Action #28 : 21/11/2025 12:34:24 p. m. : MERGE peak # 52
Action #29 : 21/11/2025 12:34:24 p. m. : SPLIT peak # 18 at 9.60

Peak Identification table :

Peak Name	RT [min]	Abs.Window[min]	Window %	Ref?	Mode
UNKNOWN_2	1.14	0.20	0.00		Nearest
heptano	1.48	0.20	0.00		Nearest
UNKNOWN_4	2.11	0.20	0.00		Nearest
UNKNOWN_5	4.10	0.20	0.00		Nearest
UNKNOWN_6	5.86	0.20	0.00		Nearest
UNKNOWN_7	6.93	0.20	0.00		Nearest
UNKNOWN_8	6.96	0.20	0.00		Nearest
UNKNOWN_9	7.71	0.20	0.00		Nearest
UNKNOWN_10	7.84	0.20	0.00		Nearest
UNKNOWN_11	8.04	0.20	0.00		Nearest
UNKNOWN_12	8.41	0.20	0.00		Nearest
UNKNOWN_13	9.32	0.20	0.00		Nearest
UNKNOWN_1	9.54	0.20	0.00		Nearest
UNKNOWN_14	10.71	0.20	0.00		Nearest

Resolve with references : No

Chromatogram : 6.2_channel1

System : 456-GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 07/11/2025 11:34:48 a. m.
Processed : 21/11/2025 12:34:23 p. m.
Printed : 21/11/2025 12:34:43 p. m.

Group Identification table :

Group Name	Group type	Parameters
------------	------------	------------

Calibration Method :

Method type : Response %

Response : Area

Standard Unit : g/l

Calibration curve file name :

Factors : Curve

Subtract internal standard mass

Response unit : Curve

Unknown mode : None

Calibration components :

Component	Level 1	Level 2	Level 3	Level 4	Control Sample
-----------	---------	---------	---------	---------	----------------

Suitability tests :

Export :

File 1 :

Export type : EXCEL

Format

Add a new sheet.

Destination :

File name : SEQNAME

File path : C:\Users\456-GC\Documents\Andres_2025\2025-11-07\MMORAN\

File ext. : .XLS

Content :

- Peak results (format :)

Post Run :

Print Method :

File name : : default_standard

Copies : 0

Summary :