


# Chromatogram : 9.1\_channel1

System : 450-GC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

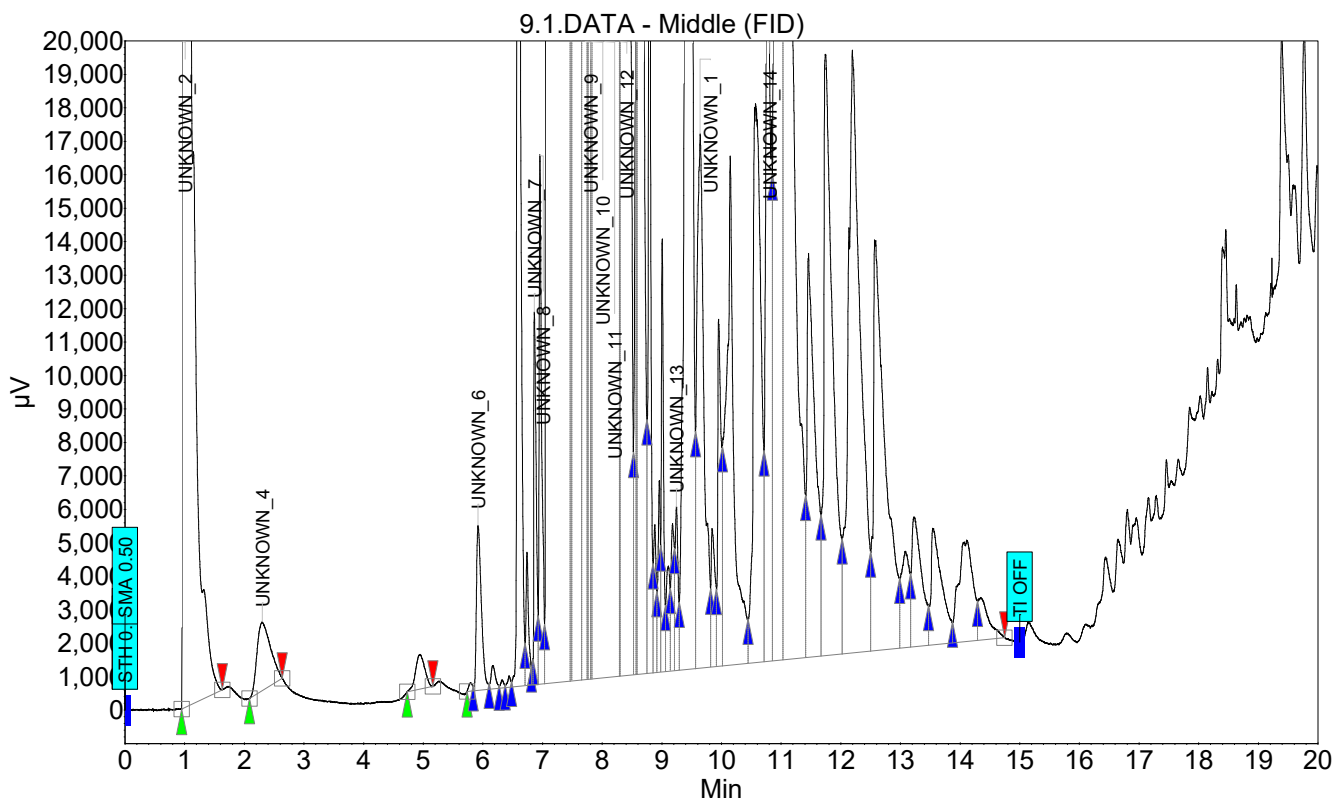

## CHROMATOGRAM METHOD REPORT :

Control method :

*Scion/Bruker/Varian 400-GC Series (Exclu*

### Autosampler (8410)

Autosampler	ENABLED
First injector used	Position 2
Syringe volume	10 µL
Advance tray	NO
Clean between injections	NO
Use injection delay	NO
Delay between injections	1.0 min
Injection mode	Std (Split/Splitless)
Sample penetration depth	90 %
Solvent penetration depth	90 %

### Abort clean

# Chromatogram : 9.1\_channel1

System : 450-CC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

Vial	I
Volume	5.0 uL
Strokes	1
Drawup speed	5.0 uL/s

Clean mode

Pre-injection solvent flushes	3
Pre-injection sample flushes	0
Post-injection solvent flushes	5
Clean solvent source	I

Injectors

Front (PWOC)

EFC 23 - Flow program	Rate (mL/min/min)	Step (mL/min)	Time (min)
	Initial	10.0	0.00
		Total time	0.00

EFC 23

Enabled OFF

Heat-only zone 1

Heater	OFF
Setpoint	50 °C

Middle (S/SL)

# Chromatogram : 9.1\_channel1

System : 450-CC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

Split event table	Time (min)	Split state	Split ratio
	Initial	ON	25

EFC 25

Enabled ON

Heat-only zone 2

Heater	ON
Setpoint	230 °C

Column oven

Stabilization time 0.50 min

Column oven zone

Heater ON

Temperature program	Rate (°C/min)	Step (°C)	Time (min)
	Initial	100	4.00
	25.0	200	8.00
	20.0	250	5.00
		Total time	23.50

Valve oven

Rear (Small valve oven)

Heat-only zone 3

Heater	OFF
Setpoint	50 °C

Columns

Front

# Chromatogram : 9.1\_channel1

System : 450-CC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

Carrier gas	Helium
Length	15.00 m
Inside diameter	200 um
Type	WCOT

## Middle

Carrier gas	Helium
Length	30.00 m
Inside diameter	320 um
Type	PLOT

Pressure mode	Constant flow
Column flow	2.20 mL/min
Pressure pulse	DISABLED

## Detectors

### Heat-only zone 4

Heater	OFF
Setpoint	50 °C

### Front (TCD)

Electronics	OFF
Time constant	Fast
Data rate	10 Hz

Filament temp. limit	390 °C
Filament temperature	50 °C

TCD event table	Time (min)	Range	Autozero	Polarity
	Initial	0.05	YES	POSITIVE

# Chromatogram : 9.1\_channel1

System : 450-CC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

## EFC 14

Enabled	OFF
Make-up (Helium) flow	25 mL/min
Reference (Helium) flow	30.0 mL/min

## Heat-only zone 5

Heater	ON
Setpoint	250 °C

## Middle (FID)

Electronics	ON
Time constant	Fast
Data rate	10 Hz

FID event table	Time (min)	Range	Autozero
	Initial	12	YES

## EFC 11

Enabled	ON
Make-up (Nitrogen) flow	28 mL/min
Combustion (H2) flow	30.0 mL/min
Combustion (Air) flow	300.0 mL/min

## Ouput ports

## Front

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Front	0

## Middle

# Chromatogram : 9.1\_channel1

System : 450-CC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

ENABLED	NO
---------	----

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Middle	0

Rear
------

ENABLED	NO
---------	----

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Rear	0

Valves
--------

Valve event table	Time (min)	Gas Sample Valve
	Initial	Fill

Miscellaneous
---------------

Start automatically when ready	NO
--------------------------------	----

Acquisition :  
Run Name : 9.1  
Run Id. : -1  
Run Time : 24.00

Integration method :  
Reduce Noise : Yes  
Spike Parameter : 1  
Use Relative Threshold ? : Yes  
Peak saturation level : min = 0   max = 0

Integration events :

Active	Time	Event	ON	Value
Yes	0.00	Set Peak Width		0.5
Yes	0.00	Set Threshold		0.05
Yes	0.00	Set Minimal Area		0.5
Yes	15.00	Turn Integration	OFF	

Manual actions :

# Chromatogram : 9.1\_channel1

System : 450 GC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

Manual integration action count : 17

Action #1 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 4 at 5.84  
Action #2 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 5 at 6.27  
Action #3 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 6 at 6.48  
Action #4 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 8 at 6.84  
Action #5 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 16 at 8.85  
Action #6 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 17 at 8.98  
Action #7 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 18 at 9.14  
Action #8 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 43 at 9.21  
Action #9 : 27/10/2025 12:44:39 p. m. : MOVE stop baseline # 4 from t=14.633 y=2368.682 to t=14.745 y=2168.203  
Action #10 : 27/10/2025 12:44:39 p. m. : MOVE stop peak # 36 from t=14.633 to t=14.747  
Action #11 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 13 at 7.80  
Action #12 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 45 at 7.83  
Action #13 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 12 at 7.75  
Action #14 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 15 at 8.57  
Action #15 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 11 at 7.49  
Action #16 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 46 at 8.29  
Action #17 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 48 at 8.58

Peak Identification table :

Peak Name	RT [min]	Abs.Window[min]	Window %	Ref?	Mode
UNKNOWN_2	1.14	0.20	0.00		Nearest
heptano	1.48	0.20	0.00		Nearest
UNKNOWN_4	2.11	0.20	0.00		Nearest
UNKNOWN_5	4.10	0.20	0.00		Nearest
UNKNOWN_6	5.86	0.20	0.00		Nearest
UNKNOWN_7	6.93	0.20	0.00		Nearest
UNKNOWN_8	6.96	0.20	0.00		Nearest
UNKNOWN_9	7.71	0.20	0.00		Nearest
UNKNOWN_10	7.84	0.20	0.00		Nearest
UNKNOWN_11	8.04	0.20	0.00		Nearest
UNKNOWN_12	8.41	0.20	0.00		Nearest
UNKNOWN_13	9.32	0.20	0.00		Nearest
UNKNOWN_1	9.54	0.20	0.00		Nearest
UNKNOWN_14	10.71	0.20	0.00		Nearest

Resolve with references : No

Group Identification table :

Group Name	Group type	Parameters
------------	------------	------------

# Chromatogram : 9.1\_channel1

System : 450-CC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

## Peak results :

Index	Name	Time [Min]	Quantity [% Area]	Height [µV]	Area [µV.Min]	Area % [%]
1	UNKNOWN 2	0.97	58.96	404397550.4	2687147.7	58.961
2	UNKNOWN 4	2.30	0.01	2051.6	532.5	0.012
3	UNKNOWN	4.94	0.00	1025.3	184.2	0.004
4	UNKNOWN	5.79	0.00	235.5	17.4	0.000
37	UNKNOWN 6	5.92	0.01	4928.2	477.2	0.010
5	UNKNOWN	6.17	0.00	678.8	60.5	0.001
38	UNKNOWN	6.32	0.00	242.5	14.5	0.000
6	UNKNOWN	6.44	0.00	338.7	20.7	0.000
39	UNKNOWN	6.58	0.06	40182.7	2690.3	0.059
7	UNKNOWN	6.74	0.00	3982.4	208.3	0.005
8	UNKNOWN	6.84	0.00	782.1	17.3	0.000
40	UNKNOWN 7	6.86	0.01	11142.4	532.5	0.012
9	UNKNOWN 8	6.96	0.02	15818.7	806.5	0.018
10	UNKNOWN	7.06	3.41	1103736.3	155216.6	3.406
11	UNKNOWN	7.48	0.02	45070.5	833.0	0.018
49	UNKNOWN	7.51	1.06	470204.5	48261.9	1.059
12	UNKNOWN 9	7.71	0.56	466407.2	25491.1	0.559
47	UNKNOWN	7.75	0.08	464631.0	3712.4	0.081
13	UNKNOWN	7.79	1.87	2732319.1	85305.6	1.872
45	UNKNOWN 10	7.82	1.50	2884024.1	68436.9	1.502
46	UNKNOWN 11	8.04	29.64	4360823.5	1350859.0	29.640
50	UNKNOWN	8.29	0.07	308123.8	3332.9	0.073
14	UNKNOWN 12	8.31	1.05	582248.3	47891.8	1.051
15	UNKNOWN	8.56	0.03	126374.6	1460.5	0.032
48	UNKNOWN	8.58	0.09	284112.9	3982.5	0.087
51	UNKNOWN	8.59	0.37	277548.2	16775.3	0.368
16	UNKNOWN	8.79	0.03	33193.7	1466.3	0.032
41	UNKNOWN	8.88	0.01	4422.4	231.5	0.005
17	UNKNOWN	8.96	0.01	5704.0	257.6	0.006
42	UNKNOWN	9.01	0.01	12923.4	545.6	0.012
18	UNKNOWN	9.11	0.00	3147.9	223.2	0.005
43	UNKNOWN	9.18	0.01	4402.1	270.6	0.006
44	UNKNOWN 13	9.24	0.01	4887.8	281.8	0.006
19	UNKNOWN	9.41	0.26	128189.3	11695.9	0.257
20	UNKNOWN 1	9.64	0.05	15962.6	2243.8	0.049
21	UNKNOWN	9.84	0.01	4135.9	321.2	0.007
22	UNKNOWN	9.95	0.02	10373.5	704.9	0.015
23	UNKNOWN	10.15	0.05	15210.4	2349.4	0.052
24	UNKNOWN	10.57	0.06	16707.8	2572.5	0.056
25	UNKNOWN 14	10.77	0.05	22586.3	2343.6	0.051
26	UNKNOWN	10.98	0.12	37694.8	5292.0	0.116
27	UNKNOWN	11.06	0.16	43693.2	7065.0	0.155
28	UNKNOWN	11.46	0.04	12071.7	1929.1	0.042
29	UNKNOWN	11.74	0.07	17965.7	3215.3	0.071
30	UNKNOWN	12.19	0.10	18013.7	4337.6	0.095
31	UNKNOWN	12.57	0.06	12274.0	2847.9	0.062
32	UNKNOWN	13.09	0.01	2848.4	448.1	0.010
33	UNKNOWN	13.23	0.02	3872.2	756.2	0.017
34	UNKNOWN	13.54	0.02	3479.8	744.3	0.016
35	UNKNOWN	14.12	0.02	3012.0	853.1	0.019
36	UNKNOWN	14.34	0.01	1267.1	261.1	0.006
Total			100.00	419022623.3	4557526.5	100.000



# Chromatogram : 9.1\_channel1

System : 450-GC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

## CHROMATOGRAM METHOD REPORT :

Control method :

*Scion/Bruker/Varian 400-GC Series (Exclu*

### Autosampler (8410)

Autosampler	ENABLED
First injector used	Position 2
Syringe volume	10 uL
Advance tray	NO
Clean between injections	NO
Use injection delay	NO
Delay between injections	1.0 min
Injection mode	Std (Split/Splitless)
Sample penetration depth	90 %
Solvent penetration depth	90 %

### Abort clean

Vial	I
Volume	5.0 uL
Strokes	1
Drawup speed	5.0 uL/s

### Clean mode

Pre-injection solvent flushes	3
Pre-injection sample flushes	0
Post-injection solvent flushes	5
Clean solvent source	I

### Injectors

### Front (PWOC)

EFC 23 - Flow program	Rate (mL/min/min)	Step (mL/min)	Time (min)
	Initial	10.0	0.00

# Chromatogram : 9.1\_channel1

System : 450 GC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

EFC 23 - Flow program	Rate (mL/min/min)	Step (mL/min)	Time (min)
		Total time	0.00

EFC 23

Enabled OFF

Heat-only zone 1

Heater	OFF
Setpoint	50 °C

Middle (S/SL)

Split event table	Time (min)	Split state	Split ratio
	Initial	ON	25

EFC 25

Enabled ON

Heat-only zone 2

Heater	ON
Setpoint	230 °C

Column oven

Stabilization time 0.50 min

Column oven zone

Heater ON

Temperature program	Rate (°C/min)	Step (°C)	Time (min)
---------------------	---------------	-----------	------------

# Chromatogram : 9.1\_channel1

System : 450-CC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

Temperature program	Rate (°C/min)	Step (°C)	Time (min)
	Initial	100	4.00
	25.0	200	8.00
	20.0	250	5.00
		Total time	23.50

Valve oven

Rear (Small valve oven)

Heat-only zone 3

Heater	OFF
Setpoint	50 °C

Columns

Front

Carrier gas	Helium
Length	15.00 m
Inside diameter	200 um
Type	WCOT

Middle

Carrier gas	Helium
Length	30.00 m
Inside diameter	320 um
Type	PLOT

Pressure mode	Constant flow
Column flow	2.20 mL/min
Pressure pulse	DISABLED

Detectors

Heat-only zone 4

# Chromatogram : 9.1\_channel1

System : 450-OC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

Heater	OFF
Setpoint	50 °C

## Front (TCD)

Electronics	OFF
Time constant	Fast
Data rate	10 Hz

Filament temp. limit	390 °C
Filament temperature	50 °C

TCD event table	Time (min)	Range	Autozero	Polarity
	Initial	0.05	YES	POSITIVE

## EFC 14

Enabled	OFF
Make-up (Helium) flow	25 mL/min
Reference (Helium) flow	30.0 mL/min

## Heat-only zone 5

Heater	ON
Setpoint	250 °C

## Middle (FID)

Electronics	ON
Time constant	Fast
Data rate	10 Hz

FID event table	Time (min)	Range	Autozero
	Initial	12	YES

# Chromatogram : 9.1\_channel1

System : 450-CC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

## EFC 11

Enabled	ON
Make-up (Nitrogen) flow	28 mL/min
Combustion (H2) flow	30.0 mL/min
Combustion (Air) flow	300.0 mL/min

## Ouput ports

### Front

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Front	0

### Middle

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Middle	0

### Rear

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Rear	0

## Valves

Valve event table	Time (min)	Gas Sample Valve
	Initial	Fill

# Chromatogram : 9.1\_channel1

System : 456-GC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

## Miscellaneous

Start automatically when ready	NO
--------------------------------	----

Acquisition :  
System : 456-GC  
Project : Aceites  
Run Name : 9.1  
Description  
Run Id. : -1  
Run Time : 24.00  
Vial : 7  
Rack : 0  
Analysis : Unknown  
Injection volume : 1.00

Scale :  
RT Min : 0.00  
RT Max : 20.00  
Y Min : -1,000.00  
Y Max : 20,000.00  
Divisor factor : 1.00  
Multiplier factor : 1.00  
Sample mass : 0.00

## Run log :

Injection report for Position 2

Injection occurred at 24/10/2025 12:29:57 p. m.

Scion/Bruker/Varian 400-GC Series (Exclu [Scion/Bruker/Varian 400-GC Series (Excluding 430)]  
Wrapper 5.0.8.31952  
Firmware 32834  
Driver 4.1.0.394

Preprocessing :  
Blank Subtract :  
File:  
N.A.

Integration method :  
Reduce Noise : Yes  
Spike Parameter : 1  
Use Relative Threshold ? : Yes  
Peak saturation level : min = 0 max = 0

## Integration events :

Active	Time	Event	ON	Value
Yes	0.00	Set Peak Width		0.5
Yes	0.00	Set Threshold		0.05
Yes	0.00	Set Minimal Area		0.5

# Chromatogram : 9.1\_channel1

System : 456-CC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

Active	Time	Event	ON	Value
Yes	15.00	Turn Integration	OFF	

Manual actions :

Manual integration action count : 17

Action #1 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 4 at 5.84  
Action #2 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 5 at 6.27  
Action #3 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 6 at 6.48  
Action #4 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 8 at 6.84  
Action #5 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 16 at 8.85  
Action #6 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 17 at 8.98  
Action #7 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 18 at 9.14  
Action #8 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 43 at 9.21  
Action #9 : 27/10/2025 12:44:39 p. m. : MOVE stop baseline # 4 from t=14.633 y=2368.682 to t=14.745 y=2168.203  
Action #10 : 27/10/2025 12:44:39 p. m. : MOVE stop peak # 36 from t=14.633 to t=14.747  
Action #11 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 13 at 7.80  
Action #12 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 45 at 7.83  
Action #13 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 12 at 7.75  
Action #14 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 15 at 8.57  
Action #15 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 11 at 7.49  
Action #16 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 46 at 8.29  
Action #17 : 27/10/2025 12:44:39 p. m. : SPLIT peak # 48 at 8.58

Peak Identification table :

Peak Name	RT [min]	Abs.Window[min]	Window %	Ref?	Mode
UNKNOWN_2	1.14	0.20	0.00		Nearest
heptano	1.48	0.20	0.00		Nearest
UNKNOWN_4	2.11	0.20	0.00		Nearest
UNKNOWN_5	4.10	0.20	0.00		Nearest
UNKNOWN_6	5.86	0.20	0.00		Nearest
UNKNOWN_7	6.93	0.20	0.00		Nearest
UNKNOWN_8	6.96	0.20	0.00		Nearest
UNKNOWN_9	7.71	0.20	0.00		Nearest
UNKNOWN_10	7.84	0.20	0.00		Nearest
UNKNOWN_11	8.04	0.20	0.00		Nearest
UNKNOWN_12	8.41	0.20	0.00		Nearest
UNKNOWN_13	9.32	0.20	0.00		Nearest
UNKNOWN_1	9.54	0.20	0.00		Nearest
UNKNOWN_14	10.71	0.20	0.00		Nearest

Resolve with references : No

Group Identification table :

Group Name	Group type	Parameters
------------	------------	------------

Calibration Method :

Method type : Response %  
Response : Area  
Standard Unit : g/l  
Calibration curve file name :

Factors : Curve  
Subtract internal standard mass  
Response unit : Curve

Unknown mode : None

# Chromatogram : 9.1\_channel1

System : 456-CC  
Method : Biodiesel FAME 2023  
User : User1

Acquired : 24/10/2025 12:29:57 p. m.  
Processed : 27/10/2025 12:44:39 p. m.  
Printed : 27/10/2025 01:33:14 p. m.

Calibration components :

Component	Level 1	Level 2	Level 3	Level 4	Control Sample
-----------	---------	---------	---------	---------	----------------

Suitability tests :

Export :  
File 1 :  
Export type : EXCEL  
Format  
Add a new sheet.  
Destination :  
File name : SEQNAME  
File path : DEFAULT  
File ext. : .XLS  
Content :  
- Peak results (format : )

Post Run :

Print Method :  
File name : : default\_standard  
Copies : 0

Summary :