


Chromatogram : 2.1_channel1

System : 450-GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

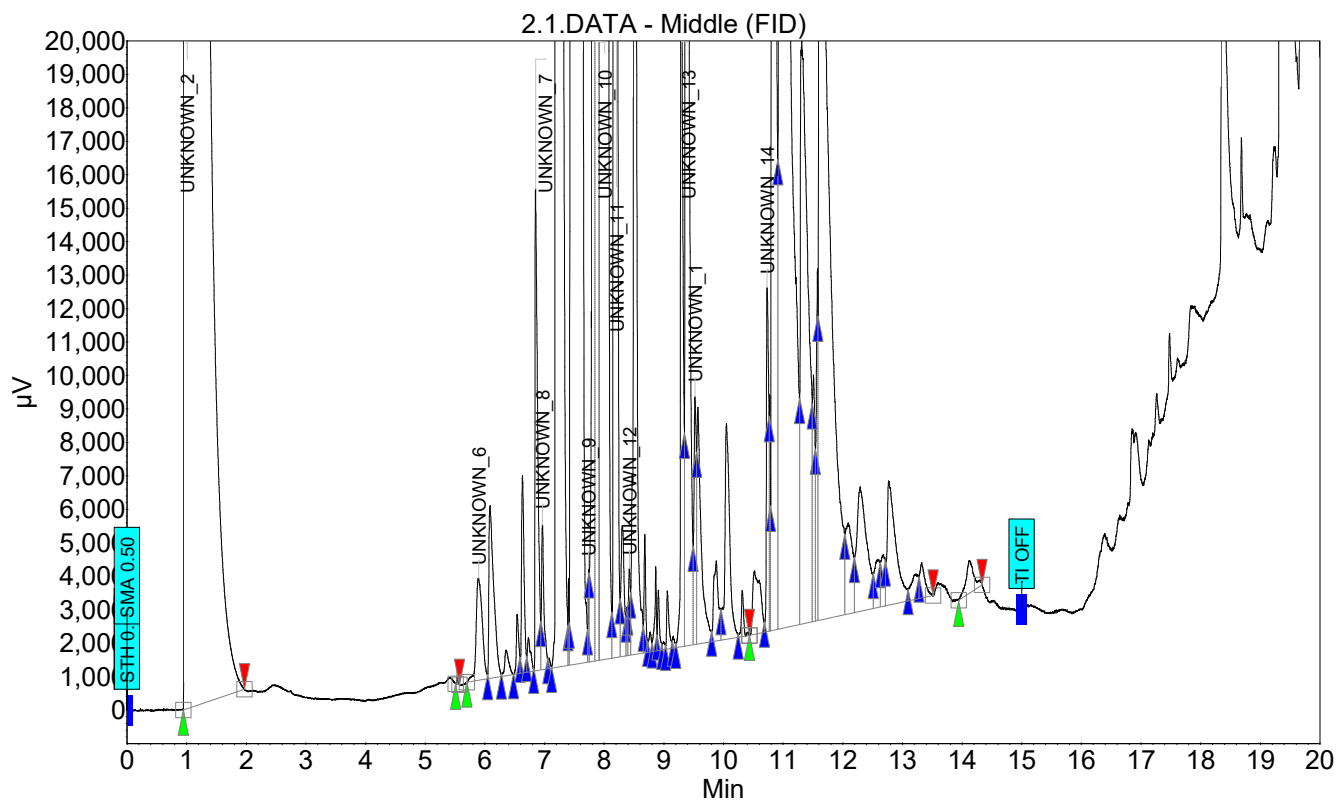

CHROMATOGRAM METHOD REPORT :

Control method :

Scion/Bruker/Varian 400-GC Series (Exclu

Autosampler (8410)

Autosampler	ENABLED
First injector used	Position 2
Syringe volume	10 uL
Advance tray	NO
Clean between injections	NO
Use injection delay	NO
Delay between injections	1.0 min
Injection mode	Std (Split/Splitless)
Sample penetration depth	90 %
Solvent penetration depth	90 %

Abort clean

Chromatogram : 2.1_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

Vial	I
Volume	5.0 uL
Strokes	1
Drawup speed	5.0 uL/s

Clean mode

Pre-injection solvent flushes	3
Pre-injection sample flushes	0
Post-injection solvent flushes	5
Clean solvent source	I

Injectors

Front (PWOC)

EFC 23 - Flow program	Rate (mL/min/min)	Step (mL/min)	Time (min)
	Initial	10.0	0.00
		Total time	0.00

EFC 23

Enabled OFF

Heat-only zone 1

Heater	OFF
Setpoint	50 °C

Middle (S/SL)

Chromatogram : 2.1_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

Split event table	Time (min)	Split state	Split ratio
	Initial	ON	25

EFC 25

Enabled ON

Heat-only zone 2

Heater	ON
Setpoint	230 °C

Column oven

Stabilization time 0.50 min

Column oven zone

Heater ON

Temperature program	Rate (°C/min)	Step (°C)	Time (min)
	Initial	100	4.00
	25.0	200	8.00
	20.0	250	5.00
		Total time	23.50

Valve oven

Rear (Small valve oven)

Heat-only zone 3

Heater	OFF
Setpoint	50 °C

Columns

Front

Chromatogram : 2.1_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

Carrier gas	Helium
Length	15.00 m
Inside diameter	200 um
Type	WCOT

Middle

Carrier gas	Helium
Length	30.00 m
Inside diameter	320 um
Type	PLOT

Pressure mode	Constant flow
Column flow	2.20 mL/min
Pressure pulse	DISABLED

Detectors

Heat-only zone 4

Heater	OFF
Setpoint	50 °C

Front (TCD)

Electronics	OFF
Time constant	Fast
Data rate	10 Hz

Filament temp. limit	390 °C
Filament temperature	50 °C

TCD event table	Time (min)	Range	Autozero	Polarity
	Initial	0.05	YES	POSITIVE

Chromatogram : 2.1_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

EFC 14

Enabled	OFF
Make-up (Helium) flow	25 mL/min
Reference (Helium) flow	30.0 mL/min

Heat-only zone 5

Heater	ON
Setpoint	250 °C

Middle (FID)

Electronics	ON
Time constant	Fast
Data rate	10 Hz

FID event table	Time (min)	Range	Autozero
	Initial	12	YES

EFC 11

Enabled	ON
Make-up (Nitrogen) flow	28 mL/min
Combustion (H2) flow	30.0 mL/min
Combustion (Air) flow	300.0 mL/min

Ouput ports

Front

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Front	0

Middle

Chromatogram : 2.1_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

ENABLED	NO
---------	----

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Middle	0

Rear

ENABLED	NO
---------	----

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Rear	0

Valves

Valve event table	Time (min)	Gas Sample Valve
	Initial	Fill

Miscellaneous

Start automatically when ready	NO
--------------------------------	----

Acquisition :
Run Name : 2.1
Run Id. : -1
Run Time : 24.00

Integration method :
Reduce Noise : Yes
Spike Parameter : 1
Use Relative Threshold ? : Yes
Peak saturation level : min = 0 max = 0

Integration events :

Active	Time	Event	ON	Value
Yes	0.00	Set Peak Width		0.5
Yes	0.00	Set Threshold		0.05
Yes	0.00	Set Minimal Area		0.5
Yes	15.00	Turn Integration	OFF	

Manual actions :

Chromatogram : 2.1_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

Manual integration action count : 30

Action #1 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 5 at 6.59
Action #2 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 30 at 6.70
Action #3 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 6 at 6.93
Action #4 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 11 at 8.27
Action #5 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 12 at 8.44
Action #6 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 13 at 8.73
Action #7 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 14 at 8.89
Action #8 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 36 at 8.97
Action #9 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 15 at 9.14
Action #10 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 19 at 10.25
Action #11 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 21 at 10.79
Action #12 : 27/10/2025 01:27:51 p. m. : MOVE valley peak # 21 from t=10.793 to t=10.790
Action #13 : 27/10/2025 01:27:51 p. m. : MOVE stop baseline # 3 from t=13.097 y=3586.222 to t=13.516 y=3425.900
Action #14 : 27/10/2025 01:27:51 p. m. : MOVE stop peak # 28 from t=13.097 to t=13.516
Action #15 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 28 at 13.09
Action #16 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 41 at 13.28
Action #17 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 27 at 12.63
Action #18 : 27/10/2025 01:27:51 p. m. : ADD peak between t=5.508 and t=5.573
Action #19 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 9 at 7.84
Action #20 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 45 at 7.91
Action #21 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 8 at 7.41
Action #22 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 12 at 8.39
Action #23 : 27/10/2025 01:27:51 p. m. : MOVE valley peak # 35 from t=8.820 to t=8.804
Action #24 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 16 at 9.34
Action #25 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 17 at 9.55
Action #26 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 21 at 10.77
Action #27 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 23 at 11.49
Action #28 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 24 at 11.58
Action #29 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 32 at 7.05
Action #30 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 9 at 7.75

Peak Identification table :

Peak Name	RT [min]	Abs.Window[min]	Window %	Ref?	Mode
UNKNOWN_2	1.14	0.20	0.00		Nearest
heptano	1.48	0.20	0.00		Nearest
UNKNOWN_4	2.11	0.20	0.00		Nearest
UNKNOWN_5	4.10	0.20	0.00		Nearest
UNKNOWN_6	5.86	0.20	0.00		Nearest
UNKNOWN_7	6.93	0.20	0.00		Nearest
UNKNOWN_8	6.96	0.20	0.00		Nearest
UNKNOWN_9	7.71	0.20	0.00		Nearest
UNKNOWN_10	7.84	0.20	0.00		Nearest
UNKNOWN_11	8.04	0.20	0.00		Nearest
UNKNOWN_12	8.41	0.20	0.00		Nearest
UNKNOWN_13	9.32	0.20	0.00		Nearest
UNKNOWN_1	9.54	0.20	0.00		Nearest
UNKNOWN_14	10.71	0.20	0.00		Nearest

Resolve with references : No

Group Identification table :

Group Name	Group type	Parameters
------------	------------	------------

Chromatogram : 2.1_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

Peak results :

Index	Name	Time [Min]	Quantity [% Area]	Height [µV]	Area [µV.Min]	Area % [%]
1	UNKNOWN 2	0.96	92.37	482960253.4	3535901.2	92.367
44	UNKNOWN	5.54	0.00	471.0	7.6	0.000
2	UNKNOWN 6	5.89	0.01	3047.2	361.5	0.009
3	UNKNOWN	6.08	0.01	5136.9	485.5	0.013
4	UNKNOWN	6.36	0.00	757.6	61.3	0.002
5	UNKNOWN	6.54	0.00	1778.3	82.7	0.002
30	UNKNOWN	6.63	0.01	5901.9	294.5	0.008
31	UNKNOWN	6.73	0.00	1012.5	64.0	0.002
6	UNKNOWN 7	6.85	0.02	14395.8	735.4	0.019
32	UNKNOWN 8	6.97	0.01	4331.5	214.9	0.006
54	UNKNOWN	7.07	0.00	308.2	12.9	0.000
7	UNKNOWN	7.19	0.48	238783.9	18534.1	0.484
8	UNKNOWN	7.40	0.00	2595.3	38.3	0.001
47	UNKNOWN	7.43	2.61	1041017.5	99981.8	2.612
9	UNKNOWN 9	7.74	0.00	2783.4	48.4	0.001
55	UNKNOWN	7.80	1.20	1037452.8	46118.4	1.205
45	UNKNOWN 10	7.85	0.90	752464.1	34387.8	0.898
46	UNKNOWN	7.92	0.02	149806.5	879.6	0.023
10	UNKNOWN	7.93	1.30	681206.6	49636.4	1.297
11	UNKNOWN 11	8.16	0.19	172236.1	7386.2	0.193
33	UNKNOWN	8.31	0.01	3897.5	242.1	0.006
12	UNKNOWN	8.38	0.00	1445.5	40.8	0.001
48	UNKNOWN 12	8.42	0.00	2591.8	99.2	0.003
34	UNKNOWN	8.49	0.06	42241.7	2464.4	0.064
13	UNKNOWN	8.68	0.00	3539.9	116.7	0.003
35	UNKNOWN	8.77	0.00	595.5	31.3	0.001
14	UNKNOWN	8.86	0.00	2524.1	82.5	0.002
36	UNKNOWN	8.91	0.00	1794.9	63.4	0.002
37	UNKNOWN	9.00	0.00	212.9	11.1	0.000
15	UNKNOWN	9.06	0.00	1731.8	76.9	0.002
38	UNKNOWN	9.16	0.00	335.4	14.5	0.000
16	UNKNOWN 13	9.29	0.03	22738.8	1067.3	0.028
49	UNKNOWN	9.37	0.09	54471.8	3330.3	0.087
17	UNKNOWN 1	9.52	0.01	7407.2	329.0	0.009
50	UNKNOWN	9.57	0.02	7062.1	605.7	0.016
18	UNKNOWN	9.88	0.01	2386.5	222.9	0.006
19	UNKNOWN	10.05	0.02	6436.6	600.2	0.016
39	UNKNOWN	10.32	0.00	1367.8	61.4	0.002
20	UNKNOWN	10.52	0.01	1891.4	277.7	0.007
21	UNKNOWN 14	10.73	0.01	10281.2	454.9	0.012
51	UNKNOWN	10.78	0.00	7058.7	133.4	0.003
40	UNKNOWN	10.85	0.07	30226.6	2866.0	0.075
22	UNKNOWN	10.96	0.23	50417.1	8747.7	0.229
23	UNKNOWN	11.31	0.07	17695.6	2593.0	0.068
52	UNKNOWN	11.51	0.01	7373.9	347.5	0.009
24	UNKNOWN	11.58	0.01	10523.0	340.7	0.009
53	UNKNOWN	11.61	0.14	27162.7	5288.0	0.138
25	UNKNOWN	12.09	0.01	2709.6	374.3	0.010
26	UNKNOWN	12.29	0.02	3721.3	730.5	0.019
27	UNKNOWN	12.58	0.00	1427.1	156.4	0.004
43	UNKNOWN	12.68	0.00	1530.9	122.4	0.003
28	UNKNOWN	12.77	0.02	3712.9	622.5	0.016
41	UNKNOWN	13.21	0.00	733.7	108.3	0.003
42	UNKNOWN	13.32	0.00	1051.7	105.3	0.003

Chromatogram : 2.1_channel1

System : 450-GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

Index	Name	Time [Min]	Quantity [% Area]	Height [μV]	Area [μV.Min]	Area % [%]
29	UNKNOWN	14.12	0.00	970.0	148.8	0.004
Total			100.00	487417009.4	3828109.5	100.000

CHROMATOGRAM METHOD REPORT :

Control method :

Scion/Bruker/Varian 400-GC Series (Exclu

Autosampler (8410)

Autosampler	ENABLED
First injector used	Position 2
Syringe volume	10 uL
Advance tray	NO
Clean between injections	NO
Use injection delay	NO
Delay between injections	1.0 min
Injection mode	Std (Split/Splitless)
Sample penetration depth	90 %
Solvent penetration depth	90 %

Abort clean

Vial	I
Volume	5.0 uL
Strokes	1
Drawup speed	5.0 uL/s

Clean mode

Pre-injection solvent flushes	3
Pre-injection sample flushes	0
Post-injection solvent flushes	5
Clean solvent source	I

Injectors

Front (PWOC)

Chromatogram : 2.1_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

EFC 23 - Flow program	Rate (mL/min/min)	Step (mL/min)	Time (min)
	Initial	10.0	0.00
		Total time	0.00

EFC 23

Enabled OFF

Heat-only zone 1

Heater	OFF
Setpoint	50 °C

Middle (S/SL)

Split event table	Time (min)	Split state	Split ratio
	Initial	ON	25

EFC 25

Enabled ON

Heat-only zone 2

Heater	ON
Setpoint	230 °C

Column oven

Stabilization time 0.50 min

Column oven zone

Chromatogram : 2.1_channel1

System : 456-OC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

Heater	ON
--------	----

Temperature program	Rate (°C/min)	Step (°C)	Time (min)
	Initial	100	4.00
	25.0	200	8.00
	20.0	250	5.00
		Total time	23.50

Valve oven

Rear (Small valve oven)

Heat-only zone 3

Heater	OFF
Setpoint	50 °C

Columns

Front

Carrier gas	Helium
Length	15.00 m
Inside diameter	200 um
Type	WCOT

Middle

Carrier gas	Helium
Length	30.00 m
Inside diameter	320 um
Type	PLOT

Pressure mode	Constant flow
Column flow	2.20 mL/min
Pressure pulse	DISABLED

Detectors

Chromatogram : 2.1_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

Heat-only zone 4

Heater	OFF
Setpoint	50 °C

Front (TCD)

Electronics	OFF
Time constant	Fast
Data rate	10 Hz

Filament temp. limit	390 °C
Filament temperature	50 °C

TCD event table	Time (min)	Range	Autozero	Polarity
	Initial	0.05	YES	POSITIVE

EFC 14

Enabled	OFF
Make-up (Helium) flow	25 mL/min
Reference (Helium) flow	30.0 mL/min

Heat-only zone 5

Heater	ON
Setpoint	250 °C

Middle (FID)

Electronics	ON
Time constant	Fast
Data rate	10 Hz

Chromatogram : 2.1_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

FID event table	Time (min)	Range	Autozero
	Initial	12	YES

EFC 11

Enabled	ON
Make-up (Nitrogen) flow	28 mL/min
Combustion (H2) flow	30.0 mL/min
Combustion (Air) flow	300.0 mL/min

Output ports

Front

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Front	0

Middle

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Middle	0

Rear

ENABLED NO

Output port program	Time (min)	Detector source	Attenuation (2^X)
	Initial	Rear	0

Valves

Valve event table	Time (min)	Gas Sample Valve
	Initial	Fill

Chromatogram : 2.1_channel1

System : 456-GC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

Miscellaneous

Start automatically when ready	NO
--------------------------------	----

Acquisition :
System : 456-GC
Project : Aceites
Run Name : 2.1
Description
Run Id. : -1
Run Time : 24.00
Vial : 1
Rack : 0
Analysis : Unknown
Injection volume : 1.00

Scale :
RT Min : 0.00
RT Max : 20.00
Y Min : -1,000.00
Y Max : 20,000.00
Divisor factor : 1.00
Multiplier factor : 1.00
Sample mass : 0.00

Run log :

Injection report for Position 2

Injection occurred at 24/10/2025 02:17:06 p. m.

Scion/Bruker/Varian 400-GC Series (Exclu [Scion/Bruker/Varian 400-GC Series (Excluding 430)]
Wrapper 5.0.8.31952
Firmware 32834
Driver 4.1.0.394

Preprocessing :
Blank Subtract :
File:
N.A.

Integration method :
Reduce Noise : Yes
Spike Parameter : 1
Use Relative Threshold ? : Yes
Peak saturation level : min = 0 max = 0

Integration events :

Active	Time	Event	ON	Value
--------	------	-------	----	-------

Chromatogram : 2.1_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

Active	Time	Event	ON	Value
Yes	0.00	Set Peak Width		0.5
Yes	0.00	Set Threshold		0.05
Yes	0.00	Set Minimal Area		0.5
Yes	15.00	Turn Integration	OFF	

Manual actions :

Manual integration action count : 30

Action #1 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 5 at 6.59
Action #2 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 30 at 6.70
Action #3 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 6 at 6.93
Action #4 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 11 at 8.27
Action #5 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 12 at 8.44
Action #6 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 13 at 8.73
Action #7 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 14 at 8.89
Action #8 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 36 at 8.97
Action #9 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 15 at 9.14
Action #10 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 19 at 10.25
Action #11 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 21 at 10.79
Action #12 : 27/10/2025 01:27:51 p. m. : MOVE valley peak # 21 from t=10.793 to t=10.790
Action #13 : 27/10/2025 01:27:51 p. m. : MOVE stop baseline # 3 from t=13.097 y=3586.222 to t=13.516 y=3425.900
Action #14 : 27/10/2025 01:27:51 p. m. : MOVE stop peak # 28 from t=13.097 to t=13.516
Action #15 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 28 at 13.09
Action #16 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 41 at 13.28
Action #17 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 27 at 12.63
Action #18 : 27/10/2025 01:27:51 p. m. : ADD peak between t=5.508 and t=5.573
Action #19 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 9 at 7.84
Action #20 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 45 at 7.91
Action #21 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 8 at 7.41
Action #22 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 12 at 8.39
Action #23 : 27/10/2025 01:27:51 p. m. : MOVE valley peak # 35 from t=8.820 to t=8.804
Action #24 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 16 at 9.34
Action #25 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 17 at 9.55
Action #26 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 21 at 10.77
Action #27 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 23 at 11.49
Action #28 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 24 at 11.58
Action #29 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 32 at 7.05
Action #30 : 27/10/2025 01:27:51 p. m. : SPLIT peak # 9 at 7.75

Peak Identification table :

Peak Name	RT [min]	Abs.Window[min]	Window %	Ref?	Mode
UNKNOWN_2	1.14	0.20	0.00		Nearest
heptano	1.48	0.20	0.00		Nearest
UNKNOWN_4	2.11	0.20	0.00		Nearest
UNKNOWN_5	4.10	0.20	0.00		Nearest
UNKNOWN_6	5.86	0.20	0.00		Nearest
UNKNOWN_7	6.93	0.20	0.00		Nearest
UNKNOWN_8	6.96	0.20	0.00		Nearest
UNKNOWN_9	7.71	0.20	0.00		Nearest
UNKNOWN_10	7.84	0.20	0.00		Nearest
UNKNOWN_11	8.04	0.20	0.00		Nearest
UNKNOWN_12	8.41	0.20	0.00		Nearest
UNKNOWN_13	9.32	0.20	0.00		Nearest
UNKNOWN_1	9.54	0.20	0.00		Nearest
UNKNOWN_14	10.71	0.20	0.00		Nearest

Resolve with references : No

Chromatogram : 2.1_channel1

System : 450-CC
Method : Biodiesel FAME 2023
User : User1

Acquired : 24/10/2025 02:17:06 p. m.
Processed : 27/10/2025 01:27:51 p. m.
Printed : 27/10/2025 01:30:58 p. m.

Group Identification table :

Group Name	Group type	Parameters
------------	------------	------------

Calibration Method :

Method type : Response %

Response : Area

Standard Unit : g/l

Calibration curve file name :

Factors : Curve

Subtract internal standard mass

Response unit : Curve

Unknown mode : None

Calibration components :

Component	Level 1	Level 2	Level 3	Level 4	Control Sample
-----------	---------	---------	---------	---------	----------------

Suitability tests :

Export :

File 1 :

Export type : EXCEL

Format

Add a new sheet.

Destination :

File name : SEQNAME

File path : DEFAULT

File ext. : .XLS

Content :

- Peak results (format :)

Post Run :

Print Method :

File name : : default_standard

Copies : 0

Summary :